





PRODUCT SPECIFICATION

| | |
|-----------------|----------------------------|
| Flow movement | Hybrid axial/radial |
| Flow direction | Down |
| Weight | 12.40 ± 0.25 g |
| Viscosity range | 1-205 cP |
| RPM range | 100-1100 RPM |
| Volume range | 500-3000 mL |
| Dimension | 57 (W) x 13 (D) x 9 (H) mm |

SAN-32

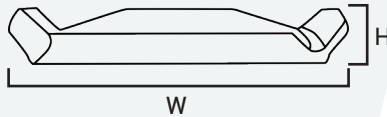

Preferred vessel :

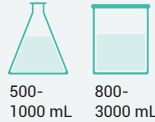

HIGH SHEAR DESIGN

Reference impeller :
Cowl disk



Utilized in emulsion and dispersion processes of :

- Emulsion and suspension polymerization
- Cosmetic creams, ceramic slurries, and polishing compounds
- Conductive inks and nanoparticle suspensions

MATERIAL TECHNICAL DATA

| | |
|-------------------------|--------------------------------|
| Shell material | polyethylene naphthalate (PEN) |
| Application temperature | -60 to 155°C |
| Color | white |
| Coefficient of friction | 0.12 |
| Hydrolysis resistance | 200 hours |
| Compliance | FDA food contact, USP <87> |
| Magnet material | SmCo permanent |
| Remanence (Br) | 9,800~11,000 Gs |
| Sterilization | gamma, steam, EtOH |

PRODUCT HIGHLIGHT



Shell made of PEN to reduce our dependency on PTFE & PFAS



Food contact safe



No cell toxicity, bio-inert



Non-hazardous materials




Estimated carbon footprint of PEN is 3.4 kg CO₂ eq/kg




Recyclable with PET bottle

FLOW DIRECTION SETTING

The up/down flow direction depends on the spin direction of the plate and the handedness of the product.

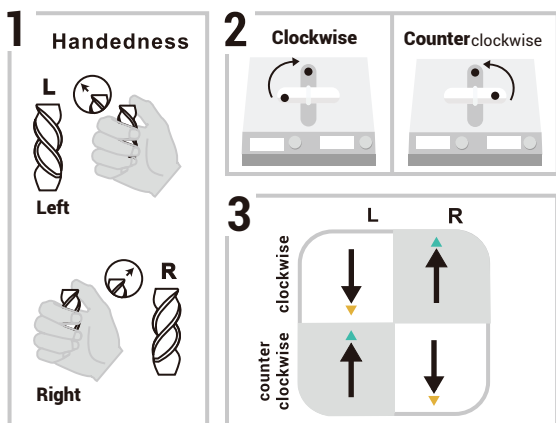

RECOMMENDED OPERATING RANGE

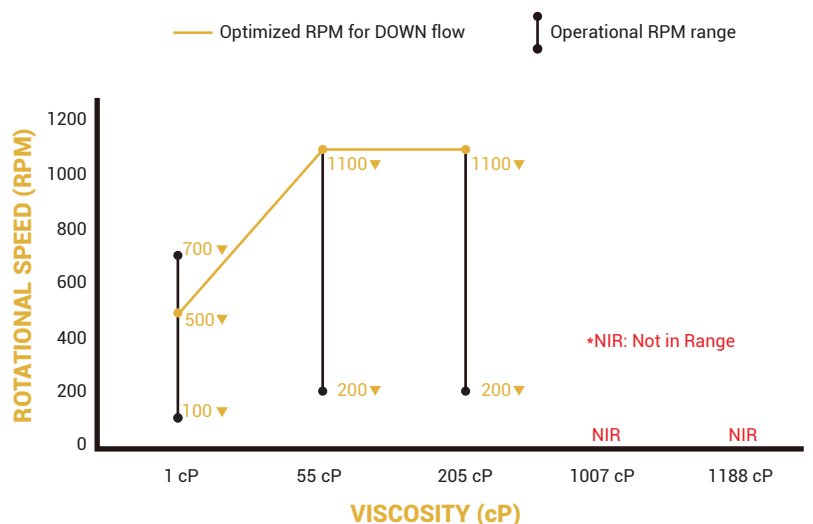

CHEMICAL RESISTANCE CHART

| | PEN | PVDF | PTFE | | PEN | PVDF | PTFE |
|---------------------------|-----|------|------|------------------------------|-----|------|------|
| Acetic acid (10%) | A | A | A | Hydrogen peroxide (30%) | A | A | A |
| Acetic acid (40%) | A | A | A | Isopropyl alcohol (IPA) | A | A | A |
| Acetone | B | D | A | Methanol | A | A | A |
| Ammonium chloride (25%) | A | A | A | Methyl ethyl ketone (MEK) | A | D | A |
| Benzene | B | A | A | Nitric acid (10%) | A | A | A |
| Chloroform | C | A | A | Nitric acid (20%) | A | A | A |
| Cyclohexane | A | A | A | Nitric acid (67%) | D | A | A |
| Cyclohexanol | A | A | A | Phosphoric acid (10%) | A | B | A |
| Calcium chloride (10%) | A | A | A | Phosphoric acid (85%) | D | B | A |
| Chromic acid (10%) | A | A | A | Potassium hydroxide (10%) | A | A | A |
| Citric acid | A | A | A | Potassium hydroxide (30%) | B | A | A |
| Diethyl ether | A | A | A | Sulfuric acid (10%) | A | A | A |
| Dimethyl sulfoxide (DMSO) | C | C | A | Sulfuric acid (30%) | C | A | A |
| Dimethyl formamide (DMF) | C | D | A | Sulfuric acid (>80%) | D | A | A |
| Ethanol | A | A | A | Sodium acetate (40%) sol. | A | A | A |
| Ethyl acetate | A | D | A | Sodium chloride (sat.) (32%) | A | A | A |
| Ethylene diamine | C | B | A | Sodium bicarbonate (sat.) | A | A | A |
| Ethylene glycol | A | A | A | Sodium hydroxide (10%) sol. | A | A | A |
| Formaldehyde (40%) | A | A | A | Sodium hydroxide (30%) sol. | A | C | A |
| Formic acid (10%) | A | A | A | Sodium hypochlorite | A | A | A |
| Formic acid (30%) | C | B | A | Styrene (monomer) | A | A | A |
| Gasoline | A | A | A | Tetrahydrofuran (THF) | B | B | A |
| Glycerol | A | A | A | Toluene | A | A | A |
| Heptane | A | A | A | Triethanolamine | C | A | A |
| Hexane | A | A | A | Vinyl chloride (monomer) | A | A | A |
| Hydrobromic acid (10%) | A | A | A | Xylene | A | A | A |
| Hydrobromic acid (47%) | A | A | A | | | | |
| Hydrochloric acid (10%) | A | A | A | | | | |
| Hydrochloric acid (37%) | A | A | A | | | | |
| Hydrofluoric acid (5%) | A | A | A | | | | |
| Hydrofluoric acid (50%) | D | B | A | | | | |

Resistance at 23°C, immersed for 30 days

Ranking definition:

A: weight change <1%, tensile > 95%


B: weight change between 1~10%, tensile > 75%

C: weight change >10%, tensile <75%

D: dissolved or swelled

SELECTION MATRIX

type of magnetic stir bars categorized by performance toward various speeds and viscosities.

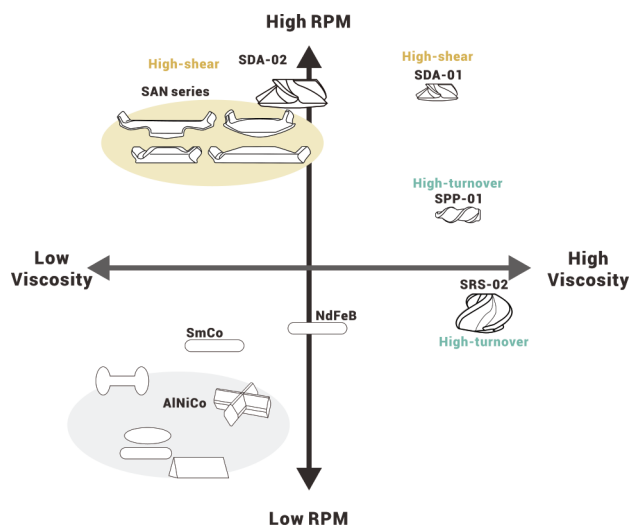

PRECAUTION

Hazard Identification:

Not classified as hazardous according to OSHA or GHS regulations.

Avoid exposure to DMSO, DMF, CHCl_3 , CCl_4 , and HF.

Check compatibility:

Use chemical resistance chart as a reference.

Conditions to Avoid:

Temperatures above 160°C.

Storage:

Magnetic stir bars should be stored in pairs to maintain their magnetic strength, away from steel or iron surfaces and other magnetic fields.

Keep them at ambient temperatures to avoid thermal stress, and handle carefully to prevent mechanical damage.

Disposal Considerations:

PEN (Polyethylene Naphthalate) can be recycled alongside PET (Polyethylene Terephthalate). Recycled as much as possible with PET or through incineration.

